


VBM Pelvic Sling

The Pelvic Sling helps to fix and stabilize fractures of the pelvis at the accident scene. Fractures of the pelvis are very painful and the bleeding resulting thereof might lead to a life-threatening condition for the patient.

The integrated inflatable cuff exerts a lateral pressure to the pelvis for best adjusting and stabilization result. This helps to reduce life-threatening inner bleeding and pain significantly.

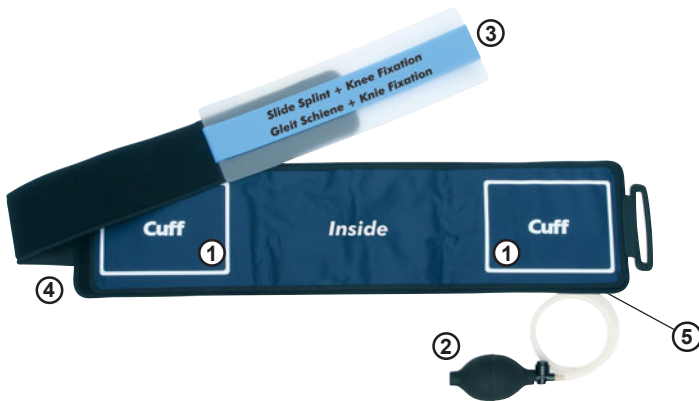

- ① **Cuffs**
Both lateral cuffs are connected and can be inflated with the hand inflator for an ideal stabilization of the pelvis.
- ② **Hand Inflator**
For inflation and deflation of both cuffs.
- ③ **Slide Splint / Knee Fixation**
To slide the sling underneath the patients pelvis. After removal used as knee fixation (see photo).
- ④ **Fixation Tape**
Length adjustable for flexible size adjustment.
- ⑤ **Stop Cock**
For a leak free transport.

Order information

for single use	
Standard, for pelvis circumference 80 - 120 cm (32" - 47")	Large, for pelvis circumference 100 - 140 cm (40" - 55")
with hand inflator	with hand inflator
REF 34-10-004	REF 34-10-005

reusable
Standard, for pelvis circumference 90 - 130 cm (36" - 51")
with manometer and hand inflator
REF 34-10-003

Manometer
0-200 mmHg
for Pelvic Sling
REF 34-18-201



VBM Medizintechnik GmbH

Einsteinstrasse 1 | DE-72172 Sulz a.N. | Germany

Tel.: +49 74 54 / 95 96 10 | Fax: +49 74 54 / 95 96 33 | e-mail: info@vbm-medical.de | www.vbm-medical.de

