

Order information



Laryngeal Tube LT, color coded syringe and biteblock

Size	Patient	Weight / Height	Color code	REF	Box
0	Newborn	< 5kg	transparent	32-01-000	1
1	Infant	5-12kg	white	32-01-001	1
2	Child	12-25kg	green	32-01-002	1
2.5	Child	125-150cm	orange	32-01-025	1
3	Adult	< 155cm	yellow	32-01-003	1
4	Adult	155-180cm	red	32-01-004	1
5	Adult	> 180cm	purple	32-01-005	1



Laryngeal Tube LT-D, color coded syringe and biteblock

Box of 10

Size	Patient	Weight / Height	Color code	REF	Box
0	Newborn	< 5kg	transparent	32-02-000-1	1
1	Infant	5-12kg	white	32-02-001-1	1
2	Child	12-25kg	green	32-02-002-1	1
2.5	Child	125-150cm	orange	32-02-025-1	1
3	Adult	< 155cm	yellow	32-02-003-1	1
4	Adult	155-180cm	red	32-02-004-1	1
5	Adult	> 180cm	purple	32-02-005-1	1

Box of 1

0	Newborn	< 5kg	transparent	32-02-100-1	1
1	Infant	5-12kg	white	32-02-101-1	1
2	Child	12-25kg	green	32-02-102-1	1
2.5	Child	125-150cm	orange	32-02-125-1	1
3	Adult	< 155cm	yellow	32-02-103-1	1
4	Adult	155-180cm	red	32-02-104-1	1
5	Adult	> 180cm	purple	32-02-105-1	1

Pediatric-Emergency Box,

LT-D #0, #1, #2, #2.5 and color coded syringe 32-02-309-1 1

Adult-Emergency Box,

LT-D #3, #4, #5, color coded syringe and biteblock 32-02-209-1 1



Laryngeal Tube LTS II, color coded syringe and biteblock

Size	Patient	Weight / Height	Color code	REF	Box
0	Newborn	< 5kg	transparent	32-05-000	1
1	Infant	5-12kg	white	32-05-001	1
2	Child	12-25kg	green	32-05-002	1
2.5	Child	125-150cm	orange	32-05-025	1
3	Adult	< 155cm	yellow	32-05-003	1
4	Adult	155-180cm	red	32-05-004	1
5	Adult	> 180cm	purple	32-05-005	1



Laryngeal Tube LTS-D, color coded syringe and one biteblock

Box of 10

Size	Patient	Weight / Height	Color code	REF	Box
3	Adult	< 155cm	yellow	32-06-003-1	1
4	Adult	155-180cm	red	32-06-004-1	1
5	Adult	> 180cm	purple	32-06-005-1	1

Box of 1

3	Adult	< 155cm	yellow	32-06-103-1	1
4	Adult	155-180cm	red	32-06-104-1	1
5	Adult	> 180cm	purple	32-06-105-1	1

Adult-Emergency Box,

LTS-D #3, #4, #5, color coded syringe and biteblock 32-06-209-1 1


VBM Laryngeal Tube

The Laryngeal Tube is a supraglottic airway device for the use in general anesthesia during spontaneous or positive pressure ventilation. In emergency it is an ideal adjunct to secure the airway during difficult airway management as an alternative technique to mask ventilation and tracheal intubation.

Securing the airway is one of the major aspects for a successful resuscitation. In the current guidelines of the European Resuscitation Council (ERC) the Laryngeal Tube has been implemented as an alternative (aid) for securing the airway during cardiopulmonary resuscitation (CPR).




Autoclavable,
Single Lumen Tube
Material: Silicone
Clean packed

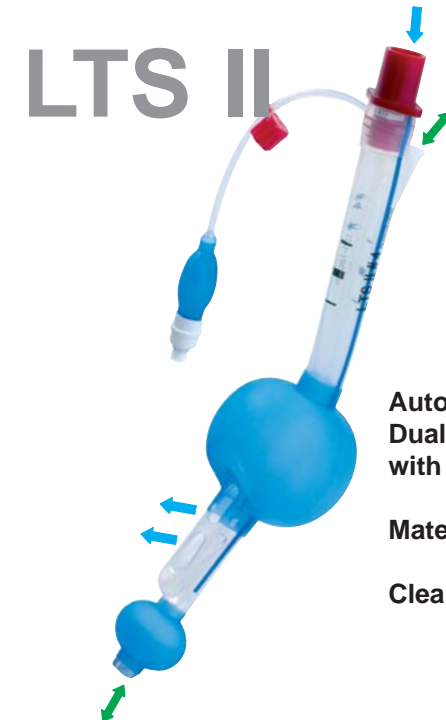

Autoclavable,
Dual Lumen Tube
with Drain Tube
Material: Silicone
Clean packed




Disposable,
Single Lumen Tube
Material: PVC
Sterile packed

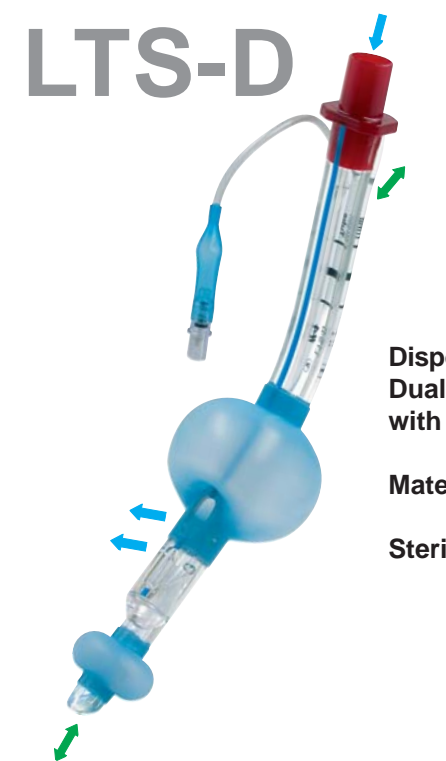


Disposable,
Dual Lumen Tube
with Drain Tube
Material: PVC
Sterile packed

Product Features

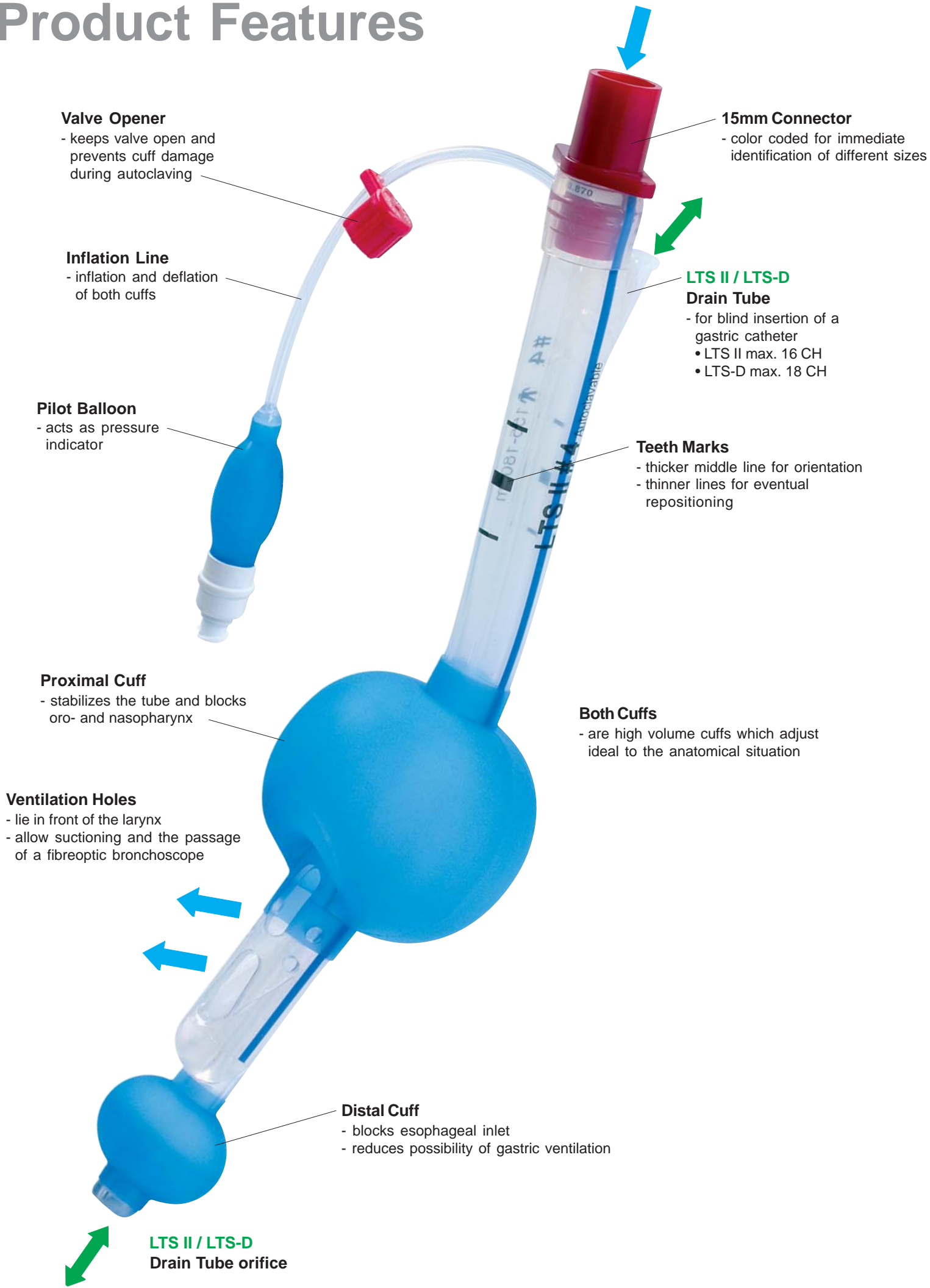

Product Features

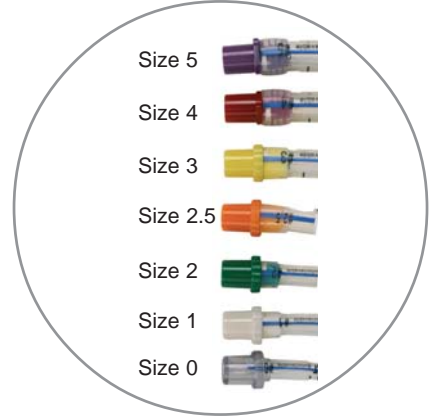


COLOR CODING

- connector and syringe are color coded
- the approximate inflation volumes are indicated on the syringe for each size to allow quick inflation of both cuffs
- in elective use a cuff pressure of 60cmH₂O is recommended (VBM Cuff Pressure Gauge "UNIVERSAL")



BITEBLOCK

- clips into the tube and protects it from occlusion
- allows tube fixation with the elastic necktape

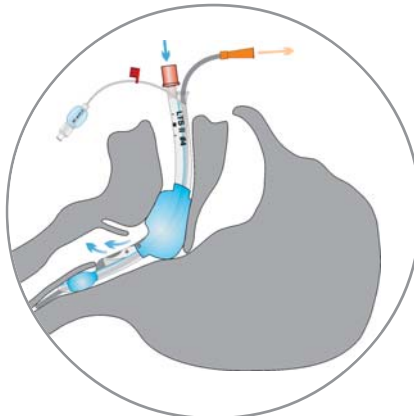

DESIGN

- form and design allow a blind insertion (without laryngoscopy) and position the Laryngeal Tube in the hypopharynx
- two Versions:
 - LT / LT-D, for ventilation
 - LTS II / LTS-D, for ventilation and drain tube

Accessories



Cuff Pressure Gauge "UNIVERSAL"

- in elective cases, to maintain the cuff pressure at 60 cmH₂O

REF 54-07-000



Syringe
with Luer connection and color code for Laryngeal Tubes

Size	LT	REF	Box
20 ml	#0-1	54-04-111	10
60 ml	#2-3	54-04-444	10
100 ml	#4-5	54-04-555	10



Biteblock with elastic necktape

- fits for Laryngeal Tube #2 and 2.5
- Size 4 **REF 40-66-040-1** for single use
- Size 5 **REF 40-76-050** reusable
- fits for Laryngeal Tube #3, 4 and 5
- Size 6 **REF 40-66-060-1** for single use
- REF 40-76-060** reusable

