

Order information



Laryngeal Tube LT, color coded syringe and biteblock

| Size | Patient | Weight/Height | Color code | REF | Box |
|------|---------|---------------|-------------|-----------|-----|
| 0 | Newborn | < 5kg | transparent | 32-01-000 | 1 |
| 1 | Infant | 5-12kg | white | 32-01-001 | 1 |
| 2 | Child | 12-25kg | green | 32-01-002 | 1 |
| 2.5 | Child | 125-150cm | orange | 32-01-025 | 1 |
| 3 | Adult | < 155cm | yellow | 32-01-003 | 1 |
| 4 | Adult | 155-180cm | red | 32-01-004 | 1 |
| 5 | Adult | > 180cm | purple | 32-01-005 | 1 |



Laryngeal Tube LT-D

| Box of 10 LT-D, color coded syringe and biteblock | REF | Box |
|---|--------------------|-----|
| 2 Child 12-25kg | green 32-02-002-1 | 1 |
| 2.5 Child 125-150cm | orange 32-02-025-1 | 1 |
| 3 Adult < 155cm | yellow 32-02-003-1 | 1 |
| 4 Adult 155-180cm | red 32-02-004-1 | 1 |
| 5 Adult > 180cm | purple 32-02-005-1 | 1 |
| Box of 1 LT-D, color coded syringe and biteblock | REF | Box |
| 2 Child 12-25kg | green 32-02-102-1 | 1 |
| 2.5 Child 125-150cm | orange 32-02-125-1 | 1 |
| 3 Adult < 155cm | yellow 32-02-103-1 | 1 |
| 4 Adult 155-180cm | red 32-02-104-1 | 1 |
| 5 Adult > 180cm | purple 32-02-105-1 | 1 |
| Emergency Box, LT-D #3, #4, #5, color coded syringe and biteblock | 32-02-209-1 | 1 |



Laryngeal Tube LTS II, color coded syringe and biteblock

| Size | Patient | Weight/Height | Color code | REF | Box |
|------|---------|---------------|-------------|-----------|-----|
| 0 | Newborn | < 5kg | transparent | 32-05-000 | 1 |
| 1 | Infant | 5-12kg | white | 32-05-001 | 1 |
| 2 | Child | 12-25kg | green | 32-05-002 | 1 |
| 2.5 | Child | 125-150cm | orange | 32-05-025 | 1 |
| 3 | Adult | < 155cm | yellow | 32-05-003 | 1 |
| 4 | Adult | 155-180cm | red | 32-05-004 | 1 |
| 5 | Adult | > 180cm | purple | 32-05-005 | 1 |



Laryngeal Tube LTS-D

| Box of 10, one color coded syringe and one biteblock | REF | Box |
|--|--------------------|-----|
| 3 Adult < 155cm | yellow 32-06-003-1 | 1 |
| 4 Adult 155-180cm | red 32-06-004-1 | 1 |
| 5 Adult > 180cm | purple 32-06-005-1 | 1 |
| Box of 1, one color coded syringe and one biteblock | REF | Box |
| 3 Adult < 155cm | yellow 32-06-103-1 | 1 |
| 4 Adult 155-180cm | red 32-06-104-1 | 1 |
| 5 Adult > 180cm | purple 32-06-105-1 | 1 |
| Emergency Box, LTS-D #3, #4, #5, color coded syringe and biteblock | 32-06-209-1 | 1 |


VBM Laryngeal Tube

The Laryngeal Tube is a supraglottic airway device for the use in general anesthesia during spontaneous or positive pressure ventilation. In emergency it is an ideal adjunct to secure the airway during difficult airway management as an alternative technique to mask ventilation and tracheal intubation.

Securing the airway is one of the major aspects for a successful resuscitation. In the current guidelines of the European Resuscitation Council (ERC) the Laryngeal Tube has been included as an additional device for airway management in cardiopulmonary resuscitation.



Autoclavable,
Single Lumen Tube
Material: Silicone
Clean packed

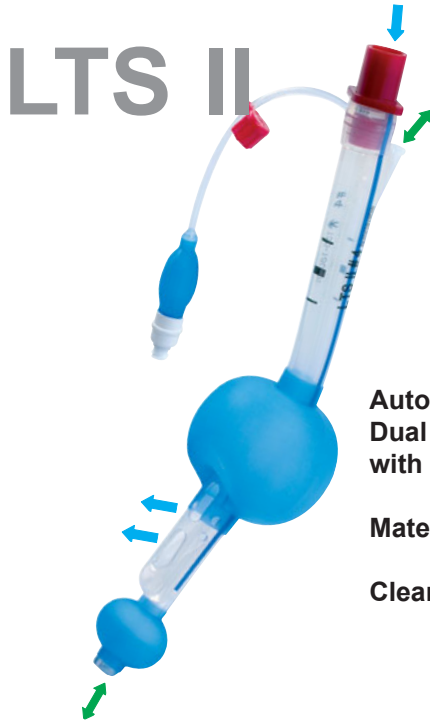

Autoclavable,
Dual Lumen Tube
with Drain Tube
Material: Silicone
Clean packed




Disposable,
Single Lumen Tube
Material: PVC
Sterile packed



Disposable,
Dual Lumen Tube
with Drain Tube
Material: PVC
Sterile packed




Product Features

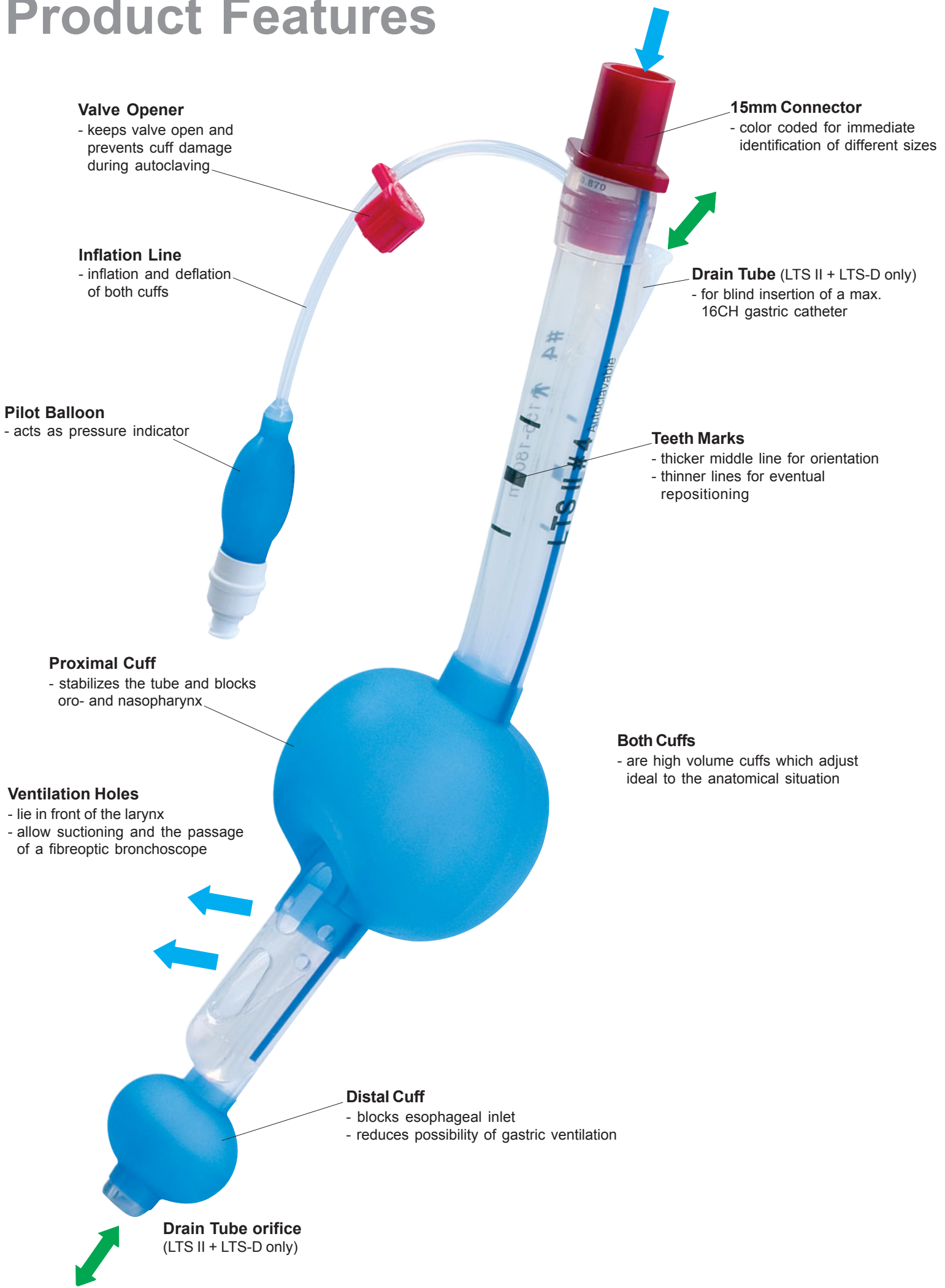

Both Cuffs
- are high volume cuffs which adjust ideal to the anatomical situation

Product Features

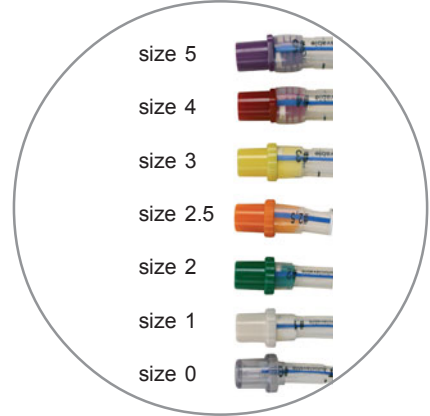


COLOR CODING

- connector and syringe are color coded
- the approximate inflation volumes are indicated on the syringe for each size to allow quick inflation of both cuffs
- in elective use a cuff pressure of 60cmH₂O is recommended (VBM Cuff Pressure Gauge "UNIVERSAL")



BITEBLOCK

- clips into the tube and protects it from occlusion
- allows tube fixation with the elastic necktape

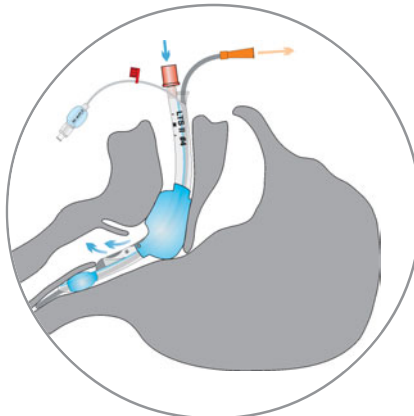

DESIGN

- form and design allow a blind insertion (without laryngoscopy) and position the Laryngeal Tube in the hypopharynx
- ease of insertion allows the Laryngeal Tube to be used by inexperienced personnel

Accessories



Cuff Pressure Gauge "UNIVERSAL"

- in elective cases, to maintain the cuff pressure at 60cmH₂O

REF 54-07-000



Syringe
with Luer connection and color code for Laryngeal Tubes

| Size | LT | REF | Box |
|-------|----------|------------------|-----|
| 20ml | size 0-1 | 54-04-111 | 10 |
| 60ml | size 2-3 | 54-04-444 | 10 |
| 100ml | size 4-5 | 54-04-555 | 10 |



Biteblock
size 6, fits Laryngeal Tube size 3,4 and 5

REF 40-70-060 autoclavable
REF 40-60-060-1 disposable

size 5, fits Laryngeal Tube size 2 and 2.5

REF 40-70-050 autoclavable
REF 40-60-050-1 disposable



Elastic Necktape
Length: 75cm

REF 40-06-075 Box 10

