

VBM Laryngobloc®

Fiber-optic Laryngoscope for single use

New concept: “One-piece“ blade and handle with reusable Light Block

■ Stability

- torsion- and break proof due to “one-piece-design”
- resistant, comparable to metal blades

■ Hygiene

- no risk of cross contamination because “one-piece-design“ prevents touch of handle with blade tip
- handle protects Light Block

■ LED light source


- extreme bright light
- constant light quality
- energy saving
- long life time

■ Cost saving

- completely disposable
- no cleaning or sterilization of the handle necessary

■ Material: Polypropylene

- excellent rigidity, no deformation
- gentle to teeth
- non-sensitive to cold temperatures

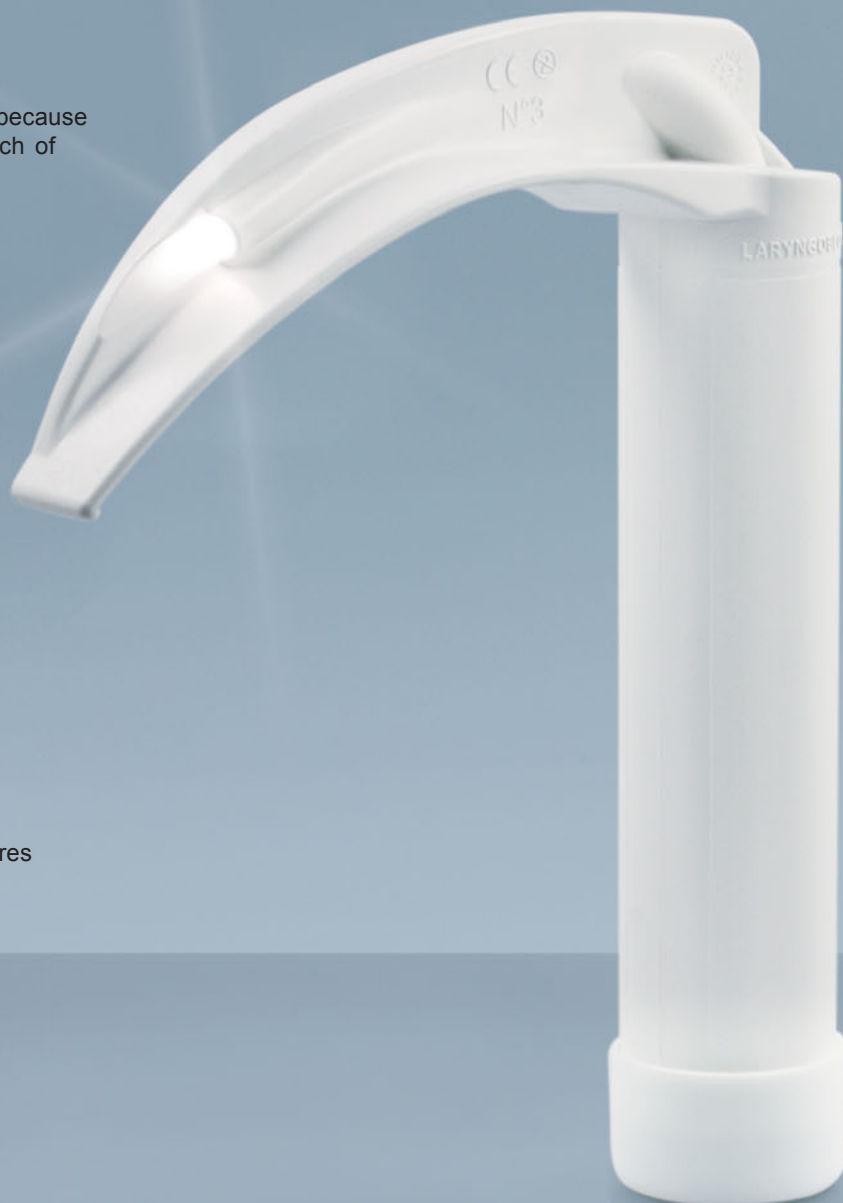

Laryngobloc®

with integrated light guide

single clean packed, disposable

Size	Description	Patient	Length	Distal Width	REF	Box
2	Macintosh	Child	95 mm	10 mm	36-20-002	20
3	Macintosh	Adult	130 mm	15 mm	36-20-003	20
4	Macintosh	Adult, large	155 mm	15 mm	36-20-004	20



Light Block

The Light Block with LED serves as light source for the Laryngobloc®.

■ Hygiene

- no cross contamination because the Light Block is protected by the handle and can be removed contamination-free after intubation

■ LED (Light Emitting Diode)

- Approximate life time: 50.000 hours

■ Battery replacement

- Typ: 2 x LR6-1.5V-AA
- Approximate life time: 15 hours in continuous use

Light Block with LED light source,
reusable, Box 1

REF 36-10-000



VBM Medizintechnik GmbH

Einsteinstrasse 1 | DE-72172 Sulz a.N. | Germany

Tel.: +49 7454 / 95 96 10 | Fax: +49 7454 / 95 96 33 | e-mail: info@vbm-medical.de | www.vbm-medical.de

